

# EquCellpen®

cordless  
mobile

Direct cold plasma therapy for  
equine medicine



# EquCellpen®

## Selected indications

### Dermatology


- Abrasion and pressure sores
- Malanders
- Mouth corner wounds
- Scars
- Saddle pressure
- Saddle bumps
- Sweet itch
- Problems with wound healing
- Burns

### Orthopaedics

- Arthritis
- Desmitis of the suspensory ligament origin (PSD)
- Desmitis of the suspensory ligament
- Kissing spines
- Tendon lesions (SDFT)
- Tendinitis


### Physiotherapy

- Trigger points
- Post-operative healing
- Haematomas
- Muscle relaxation
- Muscle spasms of the spine

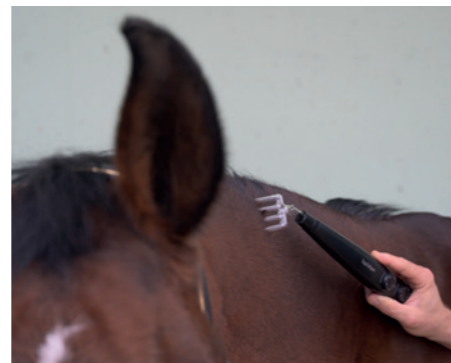


## Why the EquCellpen® is ideal for your veterinary practice

- The pen is easily manageable, easy to use and requires no auxiliary materials.
- Thanks to application-specific electrodes and variable settings, the range of its uses is extremely broad.
- With its numerous indications, the pen is suitable for day-to-day practice.
- Interruption-free work thanks to a long operating time and replaceable battery.
- And not least because it supports you in the fight against multi-resistant micro-organisms.

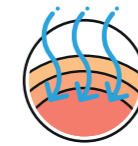

### Principle of action: direct cold atmospheric plasma (dCAP)

The cordless pen successfully generates a high voltage from the power of a small 3.7-volt battery to produce therapeutic dCAP. The plasma current can be directed precisely onto the skin surface to be treated with the help of variously shaped glass electrodes. The medical effect of dCAP has been proven in numerous studies.

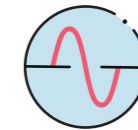

### Cell activity stimulation: promotes healing and inhibits inflammation

dCAP directly stimulates the body's cells and changes their behaviour. Healing processes are promoted through increased cell proliferation, inflammation is reduced by influencing the formation of cellular mediators, and the immune system is boosted. Tissue remodelling is improved.



### Promotion of blood circulation leads to faster healing

dCAP has an immediate and lasting hyperaemic effect on the treated tissue. The formation of new blood vessels is also stimulated. This improves metabolism and accelerates healing.



### Antimicrobial and anti-parasitic effect

The reactive species of dCAP, which contains ozone, hydrogen peroxide and other oxygen and nitrogen radicals, have a powerful antiseptic effect that does not encounter any resistance problems and is therefore effective against multi-resistant micro-organisms (e.g., MRSA). dCAP also has a fungicidal effect and has been shown to impair parasites (lice, demodex).

## The EquCellpen® and its functions



Place the interchangeable glass electrode of your choice directly on the patient's skin to transmit the electrical impulses. No need for shearing or use of a contact agent.

The Impulse/s setting determines the pulse frequency from 1 to 100 Hz according to the indication and recommendations in the Therapy Guide.

The rotary "Amplitude" control is used to switch the device on and off. Turning the control determines the device's output power and consequently the intensity of the therapy.



# EquCellpen® Set\*

Art. EQU-100



- 1 pen
- 2 glass electrodes of your choice
- 2 rechargeable batteries
- 1 charging station with cable
- 1 disinfection spray

\* The picture shows all available glass electrodes. Only two are included.

## Technical data

- Medical plasma source in accordance with DIN SPEC 91315, tested by INP (Leibniz-Institut für Plasmaforschung und Technologie e.V.)
- Cordless high-voltage generator for generating cold atmospheric plasma (CAP)
- Plasma is generated through direct dielectric barrier discharge (dDBD)
- 5 different shapes of glass electrodes with inert gas mixture
- Antibacterial device coating
- Battery insert on the back of the unit, exchangeable

Exclusive distribution partner for Switzerland and Liechtenstein

**medicalsolutiongmbh**  
the animal care company

medical solution gmbh    9500 Wil SG    Tel +41 41 741 52 75    info@medical-solution.ch  
Tonhallestrasse 47    Switzerland    Fax +41 41 741 57 75    www.medical-solution.ch

## Glass electrodes available individually in Switzerland



### Curved

**Art. CP 100-210**  
Intense, spot-focused impact with improved depth of penetration, especially for orthopaedics



### Straight

**Art. CP 100-220**  
Universal use, spot-focused or wide-area (lateral) application, dermatology and orthopaedics



### Comb

**Art. CP 100-230**  
Especially suitable for longer coats and for multi-point use: muscles, back conditions



### Disc

**Art. CP 100-240**  
Especially for dermatology, wide-area, gentler effect, for sensitive, painful areas



### Crystal

**Art. CP 100-250**  
Bundled with integrated glass crystal, for concentrated action as well as acupuncture



3-year guarantee



Made in Switzerland

Length (without electrode)	24,5 cm
Weight (without electrode)	200 g
Pulse range	1 - 100 Hz
Maximum output voltage	25 KVAC
Maximum output power	0.47 μW/cm <sup>2</sup>
Battery runtime	30 - 90 min
Battery charging time	90 min

EquCell® is a business unit of

**ActivCell Group AG**  
Luzernerstrasse 20  
6295 Mosen  
Switzerland

www.activcellgroup.com  
www.equcell.ch